
natural daylight lamp

SAD-SL01

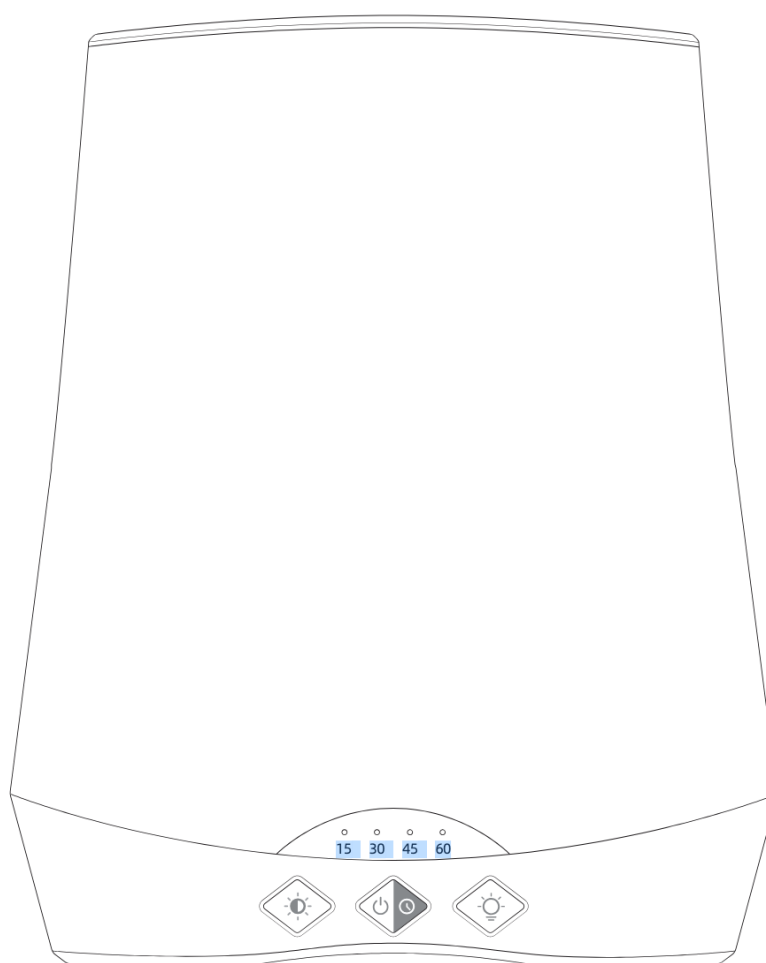Manual English

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scan voor online handleiding
scan for online manual
scanner pour le manuel française en ligne

Manual English	1
Contents	2
1. Important Safety Instructions	3
2. Features & Benefits	3
3. Pre-test	3
4. Operation	4
5. Recommended Use	5
6. Cleaning & Care	5
7. Troubleshooting	5
8. Technical Specifications	6
9. Warranty, Terms & Conditions	6

1. Important Safety Instructions

Read all instructions before use.

Warning:

- Do not use near water or in damp areas.
- Always unplug before cleaning.
- Supervise children or anyone unfamiliar with the lamp.
- Keep away from flammable materials or oxygen tanks.
- Do not use if cord or lamp is damaged.
- Indoor use only.
- Avoid dust, humidity, vibration, and direct sunlight.
- Never cover the lamp during use.
- Use only the supplied power adapter.
- Unplug during thunderstorms or long periods of non-use.

Caution:

- Certain medications or eye conditions can cause light sensitivity. Consult your doctor before use.
- Stop using immediately if you feel eye strain, headache, or discomfort.
- Do not disassemble or modify the product.
- Clean only with a soft, dry cloth. Avoid ammonia or abrasive cleaners.

2. Features & Benefits

- Natural daylight simulation (up to **10,000 LUX**)
- Adjustable **brightness** and **color temperature** (3000K–6500K)
- Timer settings: **15, 30, 45, or 60 minutes**
- Compact, portable design
- Helps improve mood, focus, and energy levels

3. Pre-test




To help determine if you have a sensitivity to strong light, we recommend that you perform this pre-test before the first use of this product.

- Place the lamp about **18 inches / 45 cm's** from your face (do not look directly into it).
- Turn it on for **10 minutes**.
- If you feel discomfort, stop using and consult a doctor.
- If no discomfort, you can gradually increase session time or move the lamp closer.

4. Operation


- Plug the adapter into the lamp and a wall outlet.
- **Power On/Off:** Hold the timer button for 1.5 seconds.
- Position the lamp as close as 6 inches/15 cm, but within 24 inches/60 cm of your face and should be slightly off-center (see diagram).
- **Adjust Brightness**  : Tap or hold the Light Intensity switch to change between 6 brightness levels or stepless mode (960–10,000 LUX).
- **Adjust Color Temperature**  : Hold the Color Temperature switch to change between **3000K (warm)** and **6500K (cool)**.
- **Set Timer**  : Press to select **15, 30, 45, or 60 minutes**. If no timer is set, the light stays on.
- Position the lamp **6–24 inches / 15–60 cm** from your face and slightly off-center. You may read, work, or relax during sessions.

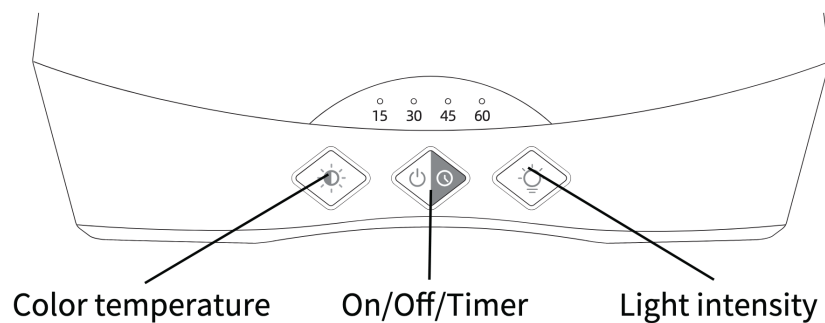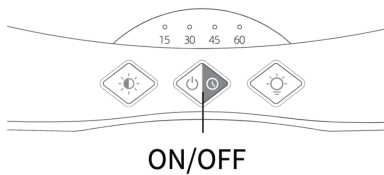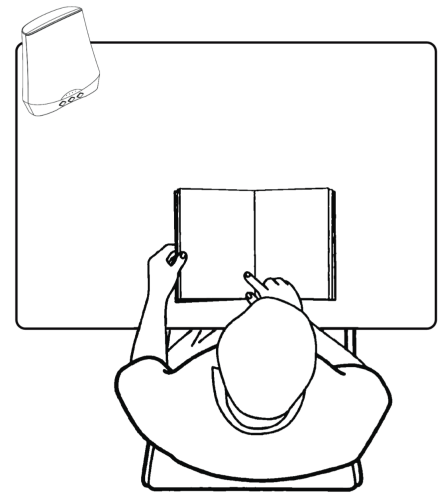5. Recommended Use

- Use **30–60 minutes daily**, preferably in the morning.
- Consistent daily use provides the best results.
- If you feel overstimulated, reduce light intensity or session length.
- If you feel sluggish, increase light intensity or sit closer.
- Sessions can be split into shorter intervals throughout the day.

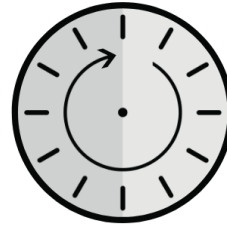30 Minutes
to
1 hour/day

6. Cleaning & Care

- Turn off, unplug, and let the lamp cool before cleaning.
- Wipe with a **soft, dry cloth**.
- Do **not** use liquids, solvents, or abrasive materials.
- Ensure the lamp is completely dry before plugging it back in.

7. Troubleshooting

- Ensure the power cord is fully connected.
- Try another wall outlet.
- Do not use with dimmers, timers, or extension cords.

8. Technical Specifications

Item	Specification
Model	KHYR SAD-SL01
Size	203 × 164 × 41.3 mm
Power	12W
Adapter	AC 110–240V → DC 12V/2A
Light Intensity	960–10,000 LUX
Color Temperature	3000K–6500K
Timer	15 / 30 / 45 / 60 min

9. Warranty, Terms & Conditions

The KHYR SAD-SL01 comes with a **12-month warranty** covering defects in materials and workmanship.

To make a warranty claim go to www.khyr.eu/garantie and fill in the form. For questions regarding warranty go to www.khyr.eu/contact or email to support@khyr.eu.

The warranty becomes void in the following situations:

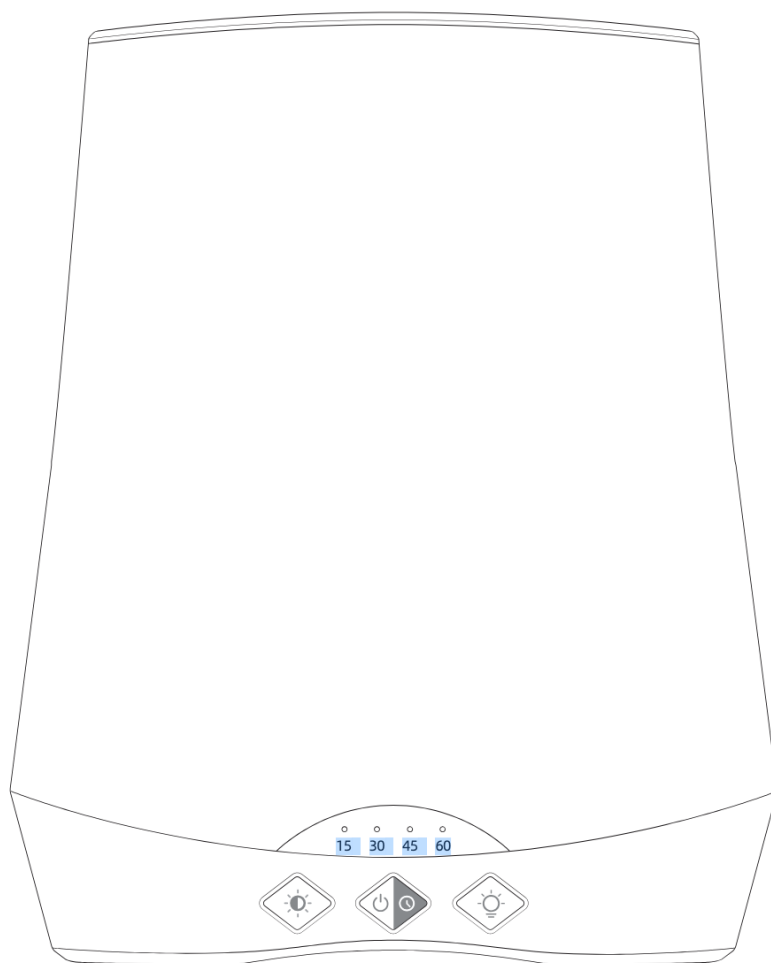
- Incorrect installation or use as described in this manual
- Accidental or intentional damage of the device
- Use with incompatible products or voltages
- Operation outside stated specifications
- Opening the Repairs or modifications by unauthorised persons

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