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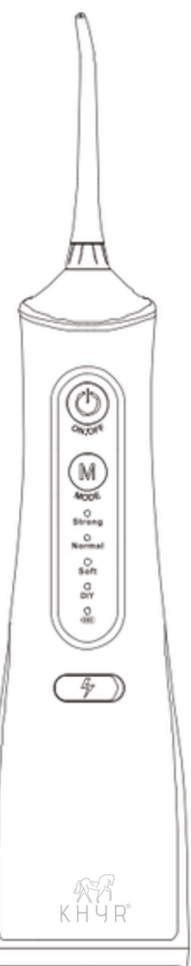
Before using the product, please read this manual carefully and keep it for future reference.

Safety Instructions


Attention

- If you do not understand the operating procedures and safety instructions, do not use this product to avoid injury caused by improper operation.
- Do not operate the product by any means other than those described in this manual.
- Do not use the water flosser product when it is plugged into a power source.
- Make sure the USB-C port is properly sealed with the rubber sealer at all times, except during charging, to prevent water and dirt from entering.
- Do not disassemble this product or modify its components to avoid damage, high-voltage electric shock, or invalidating the warranty.
- Do not use water with a temperature above 40°C in the water tank or for cleaning the mouth to avoid burns.
- If the water flosser falls into water while charging, unplug the power cord before retrieving the product.
- If the power cord or plug is damaged, replace it before use.
- Do not expose the product to fire or heat, and do not charge, use, or store it in high temperatures.

- Do not use this product for purposes other than oral cleaning.
- Place this product in a safe dry location.
- Do not immerse this product in water to prevent malfunction, fire, electric shock, explosion, or other hazards.
- Before connecting the power plug, check whether the voltage marked on the power adapter matches the local voltage.
- Do not damage or modify the power cord. Do not stretch, bend, or twist it forcefully, place heavy objects on it, or trap it between objects.
- Do not place the product, especially the power cord, near hot objects.
- Unplug the power plug and charging plug before performing maintenance.
- Do not store the product in high temperatures or humid environments, and avoid direct sunlight for longer time.
- If not used for a long time, unplug the power plug and charging plug to prevent fire or danger.

- Some mouthwash products may corrode this product; do not use mouthwash in this device.
- This product is not suitable for children, infants, or individuals who cannot operate it without supervision or instructions from a responsible person. Keep it out of the reach of children and babies. Do not use it as a toy.
- Patients with periodontal disease or those who have undergone oral surgery in the past two months should consult a doctor before use.
- When plugging in or unplugging the power supply, always hold the plug instead of the power cord.

Function overview

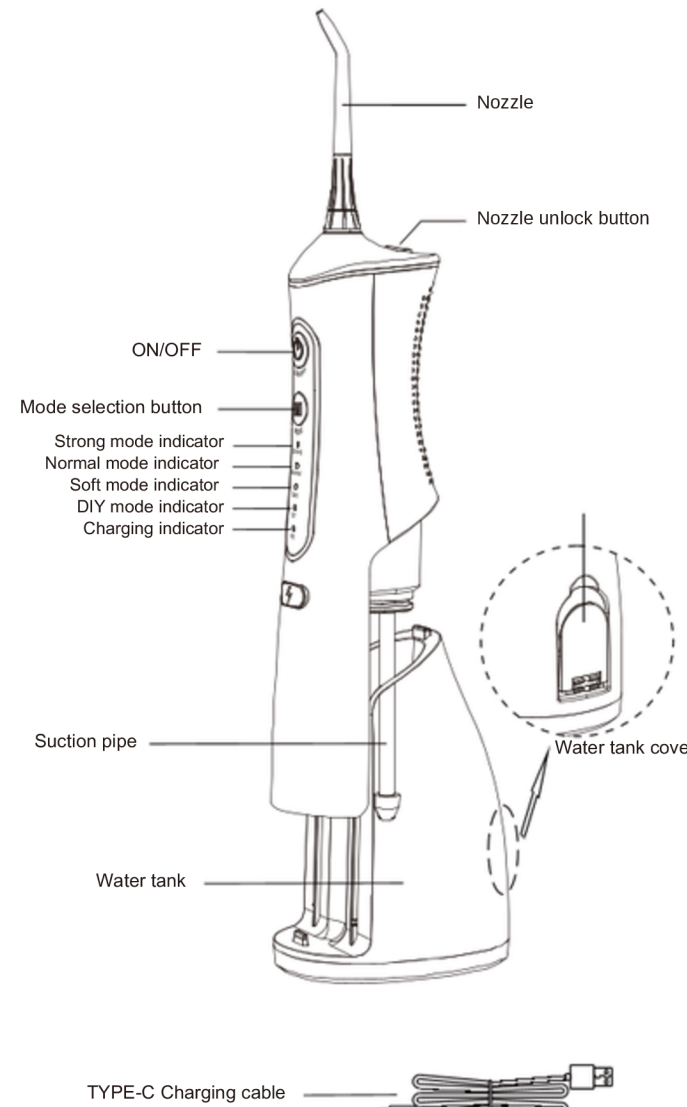Operation instructions

Before use, please charge the product as follows:

- Ensure that both the product and your hands are dry before connecting the USB plug to a 5V adapter (provided by the user) and plugging it into the socket (Figure 1).
- Remove the charging cover, insert the other end of the USB cable into the charging port of the device, and ensure the charging plug is properly inserted (Figure 2).

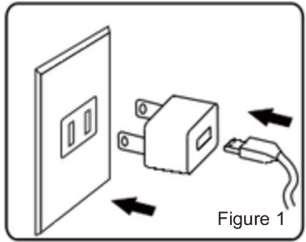
Figure 1

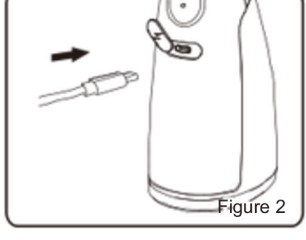
Figure 2

TIP
If the indicator light flashes during or after use, the battery is low and needs to be charged. It takes 2-3 hours to fully charge. While charging, the indicator light flashes to indicate the charging process. The indicator light stops flashing when the battery is fully charged.

When using, please follow the steps below:
Installing or replacing the nozzle: Press and hold the button on top of the main unit to unlock the nozzle slot. Then, insert the nozzle into the slot, ensuring it is properly in place (Figure 3-1). To replace the nozzle, press and hold the button again while pulling the nozzle out (Figure 3-2).

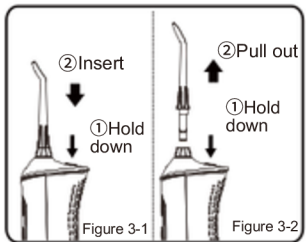
Figure 3-1


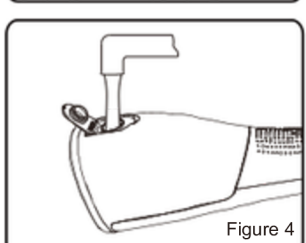
Figure 3-2


- Fill the water tank with water.
- There are two methods for filling the water tank:
 - Open the water tank cover, hold the water flosser horizontally, fill the tank with water, and then close the cover (Figure 4).
 - Pull down the water tank, remove it and hold it vertically, fill it with water, then push it back up to install (Figure 5).

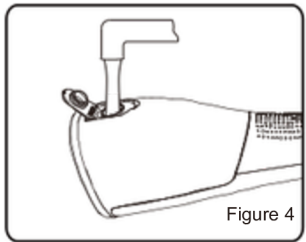
Figuur 4

The water tank must be pushed into place without leaving gaps to prevent leaks.


- Select the desired mode by pressing the mode button. You can also switch to a different mode by pressing the mode button during use.

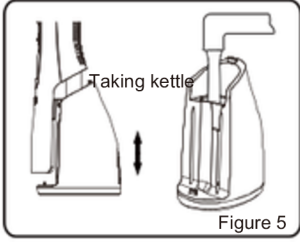
Figuur 5

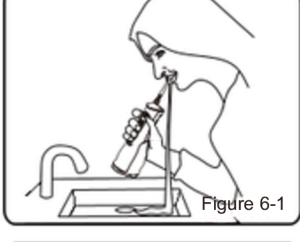Figuur 6-1

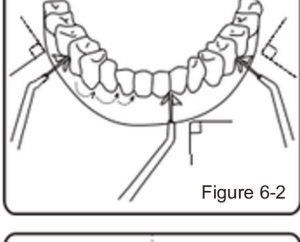Figuur 6-2

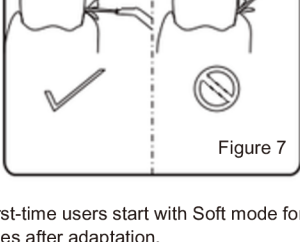Figuur 7

- 4 water pressure modes**
This product has four modes: **Strong**, **Normal**, **Soft**, and **DIY** mode. Their specific functions are as follows:
- Strong mode: Deep and fast cleaning.
 - Normal mode: Standard water pressure for daily cleaning.
 - Soft mode: Gentle cleaning, suitable for people with sensitive teeth.
 - DIY mode: A custom setting mode that allows you to freely adjust the water pressure.

DIY Mode Usage:

- Switch to DIY mode and press and hold the ON/OFF button. The water pressure will cycle from weak to strong.
- Pay attention to the change in water pressure. When you feel that the pressure reaches a comfortable or suitable level, release the button.
- The selected water pressure is now set for DIY mode, and you can start using the device.
- The product has a memory function that saves your DIY water pressure setting for the next use.

TIP:

Dental professionals recommend that first-time users start with Soft mode for one week before switching to other modes after adaptation.
1. Before use, hold the water flosser vertically, ensure the nozzle extends into the oral cavity, and align it with the teeth or gums. Keep your mouth slightly open to allow the water to flow out smoothly (Figure 6-1 & 6-2).

Routine maintenance

To prevent the spread of bacteria, ensure that no residual water remains in the dental flosser. If you do not plan to use it for a long time, clean the product and store it after drying.

Use clean water or a neutral detergent to clean the product. Do not use corrosive or abrasive cleaners, as they may damage the device. Do not rinse with hot water above 50°C.

- Clean the body with a dry cloth or a damp cloth.
- Do not immerse the device in water for cleaning.

Cleaning the water tank:

- Clean it with a dry or damp cloth, or remove it and rinse with clean water.
- If you do not plan to use it for a long time, make sure to empty and clean the water inside.

Cleaning the nozzle:

- Clean it with a dry or damp cloth, or remove it and rinse with clean water. Then, wipe it dry with a cloth.
- Do not bend, pull, or twist the nozzle hose.
- Dental experts recommend replacing the nozzle every 6 months.

Troubleshooting

Problem	Possible cause	Solution
Device stops functioning	New purchase or long time no charge leads to low battery	Charge for 120mins
Low water pressure during use or no water spray	Something is blocked in the nozzle	Clean or replace the nozzle
	Suction inlet is blocked	Clean the suction inlet
Only works for a short time after charging	Device has been idle for a while	Pull out the nozzle and flush the nozzle hole
	The battery has reached the end of its useful life	Contact your supplier for replacement
The charging indicator does not light up during charging	Insufficient charging time	Charge for more than 2hrs
	Power connector not connected properly	Connect power supply properly
	Product charging circuit is faulty	Contact your supplier

Waste disposal

This product uses rechargeable batteries. At the end of its lifespan, before discarding the product, be sure to remove the battery and recycle or dispose of it in accordance with local waste management regulations. If in doubt, please consult your local waste management department.

Support Contact
MCIYSS
Laan van Waalhaven 139A
2497 GK Den Haag
support@khyr.eu
www.khyr.eu